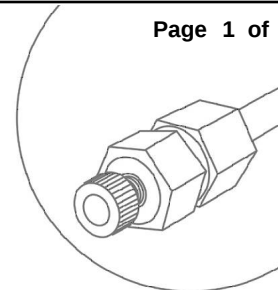



Cation-Exchange Std. (AL0-3040)

Column: Luna 5u SCX 100A**Dimensions:** 250 x 4.6 mm ID**Order No:** 00G-4398-E0**Elution Type:** Isocratic**Eluent A:** 0.15M (NH₄)₂H₂PO₄, pH 6.0**Pct. A:** 100**Flow Rate:** 1 mL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs. - Variable Wave. @ 254 nm (ambient)

Page 1 of 1

**Products used in this Application:****LUNA[®]**

Analyst Note: This standard is used to test the following Phenomenex silica-based ion-exchange columns:
SCX, SA, CM.
NOTE: Separation conditions and results may vary for your column. Contact your Technical Representative.

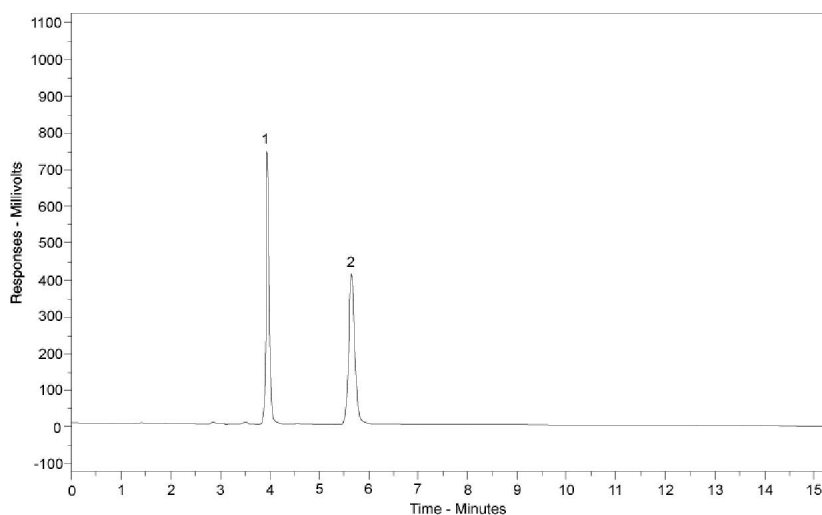**ANALYTES:**

- 1 Uracil
- 2 Cytosine

www.phenomenex.comPhenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.