



Solvents (mix 1)

Column: Zebron ZB-5MS
Phase: 5% Polysilarylene - 95% Polydimethylsiloxane copolymer
Dimensions: 30 meters x 0.25 mm x 0.25 µm
Order No: 7HG-G010-11
Oven Profile: 40°C for 2 min to 70°C at 4°C/min to 280°C at 25°C/min
Carrier Gas: Constant Flow Helium, 1.1 mL/min
Injection: split 100:1 0.2 µL @ 250 °C
Detection: Mass Selective (MSD) (240 °C)

Page 1 of 2

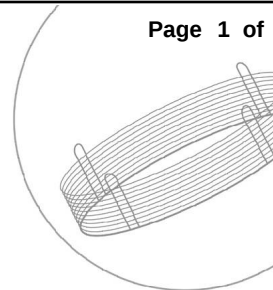

Products used in this Application:

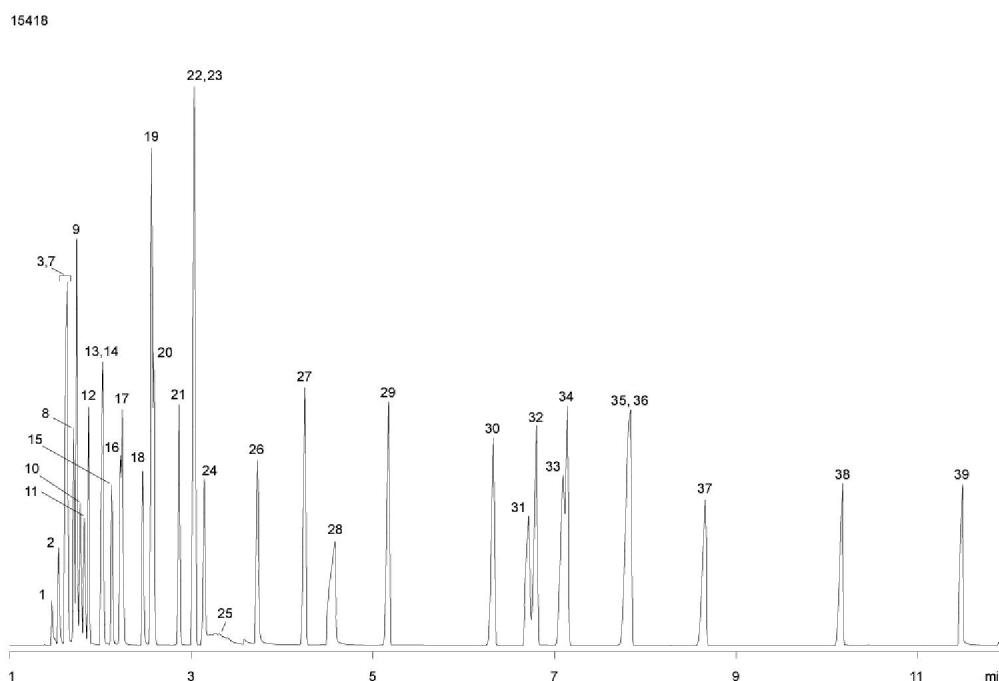ANALYTES:

- | | | |
|----------------------|--------------------------|----------------------|
| 1 Methanol | 10 Allyl Alcohol | 19 Benzene |
| 2 Ethanol | 11 n-Propanol | 20 Cyclohexane |
| 3 Acetonitrile | 12 MTBE | 21 iso-Octane |
| 4 Acetone | 13 2-Butanone | 22 Heptane |
| 5 Pentane | 14 Hexane | 23 Trichloroethylene |
| 6 Isopropanol | 15 Ethyl Acetate | 24 1,4-Dioxane |
| 7 Ethyl Ether | 16 Chloroform | 25 Ethylene glycol |
| 8 tert-Butanol | 17 Tetrahydrofuran (THF) | 26 Pyridine |
| 9 Methylene chloride | 18 1,2-Dimethoxy ethane | 27 Toluene |

©Copyright 2000-2007 Phenomenex Inc.

For more information call your Phenomenex Representative at (310) 212-0555



www.phenomenex.com

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.


 ...breaking with tradition™

USA
 tel.: (310) 212-0555
 fax: (310) 328-7768
 email: info@phenomenex.com

Puerto Rico
 (800) 541-HPLC
 (310) 328-7768
 info@phenomenex.com

Canada
 (800) 543-3681
 (310) 328-7768
 info@phenomenex.com

France
 01 30 09 21 10
 01 30 09 21 11
 franceinfo@phenomenex.com

United Kingdom
 01625-501367
 01625-501796
 ukinfo@phenomenex.com

Germany
 06021-58830-0
 06021-58830-11
 anfrage@phenomenex.com

New Zealand
 09-4780951
 09-4780952
 info@phenomenex.co.nz

Australia
 02-9428-6444
 02-9428-6445
 info@phenomenex.com.au

Solvents (mix 1)

Page 2 of 2

- 28 N,N-Dimethylformamide
- 29 Tetrachloroethylene
- 30 Chlorobenzene
- 31 3-Methylpyridine
- 32 Ethylbenzene
- 33 m-Xylene
- 34 p-Xylene
- 35 Styrene
- 36 o-Xylene
- 37 Anisole
- 38 Benzaldehyde
- 39 N-Methylpyrrolidinone

www.phenomenex.comPhenomenex products are available worldwide. For the distributor in your country,
contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.

USA
tel.: (310) 212-0555
fax: (310) 328-7768
email: info@phenomenex.com

Puerto Rico
(800) 541-HPLC
(310) 328-7768
info@phenomenex.com

Canada
(800) 543-3681
(310) 328-7768
info@phenomenex.com

France
01 30 09 21 10
01 30 09 21 11
franceinfo@phenomenex.com

United Kingdom
01625-501367
01625-501796
ukinfo@phenomenex.com

Germany
06021-58830-0
06021-58830-11
anfrage@phenomenex.com

New Zealand
09-4780951
09-4780952
info@phenomenex.co.nz

Australia
02-9428-6444
02-9428-6445
info@phenomenex.com.au