


Food Additives

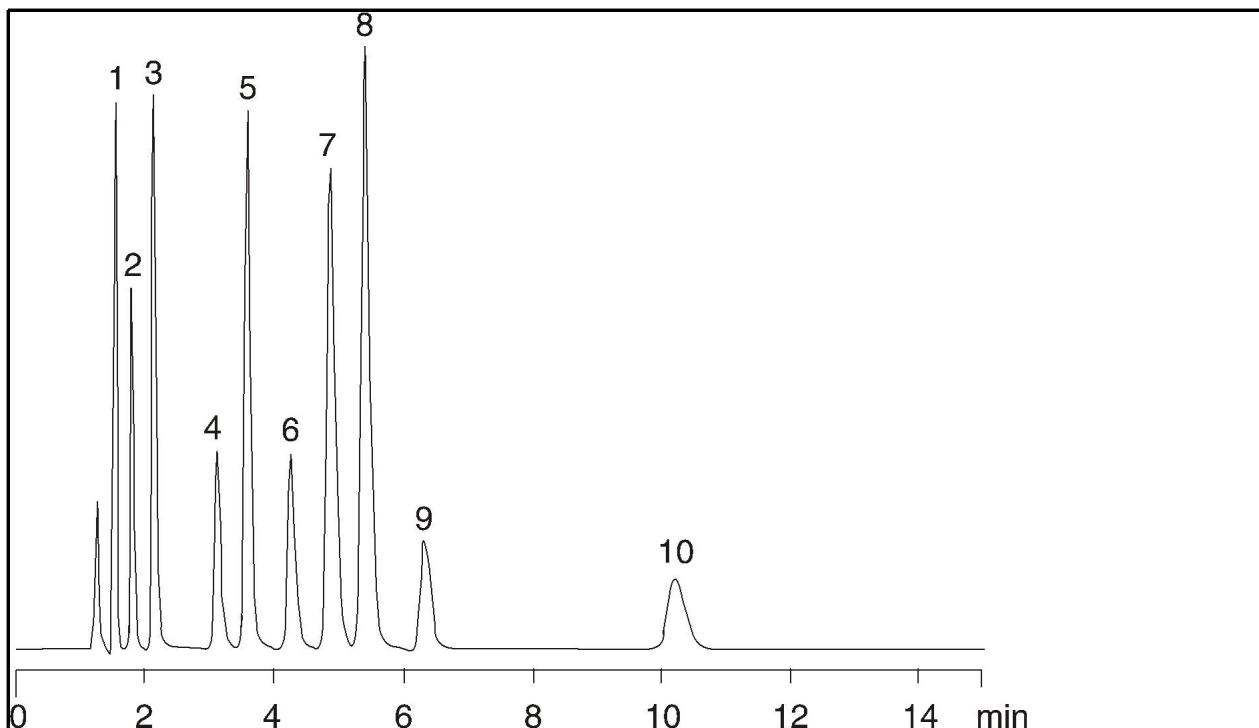

Column: Prodigy 5 μ m C6H5 / Phenyl PH-3

Dimensions: 150 x 4.6 mm ID

Order No:

Elution Type: Isocratic

Eluent A: Acetonitrile

Pct. A: 25

Eluent B: 20 mM Potassium phosphate pH 3.0

Pct. B: 75

Flow Rate: 1.5 mL/min

Temperature: 40 °C

Detection: UV-Vis Abs. - Variable Wave. @ 230 nm / 210 nm

SAMPLE ANALYTE

- 1 Saccharin
- 2 Aspartame
- 3 *p*-Hydroxybenzoic acid
- 4 Caffeine
- 5 Benzoic acid
- 6 *p*-Hydroxybenzoic acid methyl ester
- 7 *p*-Toluic acid
- 8 Dehydroacetic acid
- 9 *p*-Hydroxybenzoic acid ethyl ester
- 10 *p*-Hydroxybenzoic acid propyl ester

SPE:

For more information call your Phenomenex Representative at (310) 212-0555

Copyright 2000 Phenomenex Inc.

Ref. No.: 10257/BRO/5