

Soft Drink Additives

Column: Synergi 4u Polar-RP 80A

Dimensions: 250 x 4.6 mm ID

Order No: 00G-4336-E0

Elution Type: Isocratic

Eluent A: 20 mM Potassium Phosphate, pH 2.3

Pct. A: 55


Eluent B: Methanol

Pct. B: 45

Flow Rate: 1 mL/min

Col. Temp.: ambient


Detection: UV-Vis Abs. - Variable Wave. @ 220 nm (ambient)

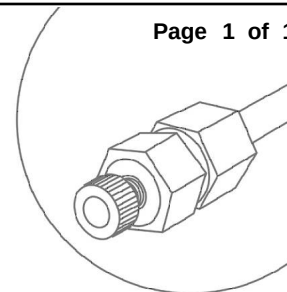

Page 1 of 1

Products used in this Application:

synergi™

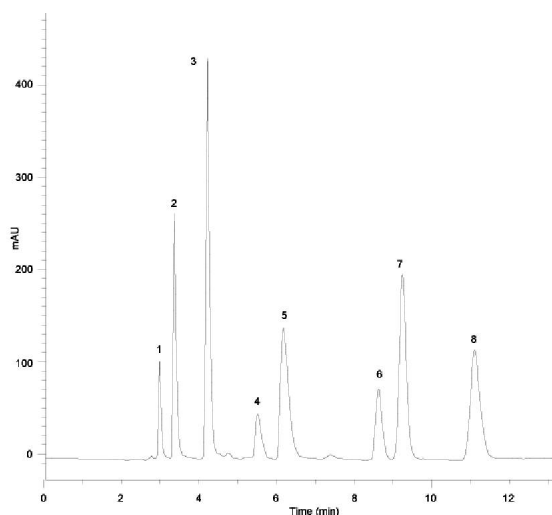

ANALYTES:

- 1 Ascorbic Acid
- 2 Acesulphame K
- 3 Saccharin
- 4 Aspartame
- 5 Quinine
- 6 Sorbate
- 7 Benzoate
- 8 Caffeine



www.phenomenex.com

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.