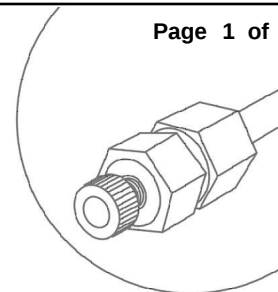
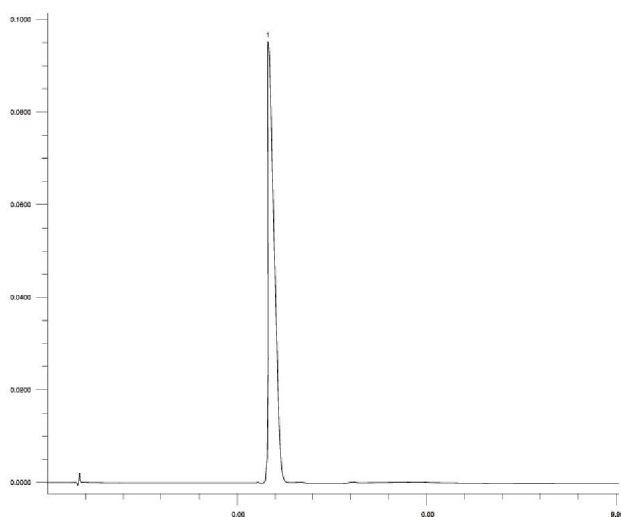


Coenzyme Q10

Column: Luna 5u C18(2) 100A**Dimensions:** 150 x 4.6 mm ID**Order No:** 00F-4252-E0**Elution Type:** Isocratic**Eluent A:** Ethanol**Pct. A:** 65**Eluent B:** Methanol**Pct. B:** 35**Flow Rate:** 1 mL/min**Col. Temp.:** 25 °C**Detection:** UV-Vis Abs. - Variable Wave. @ 275 nm (ambient)

Page 1 of 1

**Products used in this Application:****LUNA®****ANALYTES:**

- 1 Coenzyme Q10

**www.phenomenex.com**

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.