



Butter Triglycerides on ZB-5HT 15m x 0.32mm x 0.10µm

Page 1 of 1

Column: Zebron ZB-5HT
Phase: 5%-diphenyl-95%-dimethyl polysiloxane copolymer
Dimensions: 15 meters x 0.32 mm x 0.1 µm
Order No: 7EM-G015-02
Oven Profile: 100°C to 400°C at 14°C/min for 10 min
Carrier Gas: Constant Flow Helium, 1.8 mL/min
Injection: On-Column 2 µL @ 103 °C
Detection: Flame Ionization (FID) (400 °C)

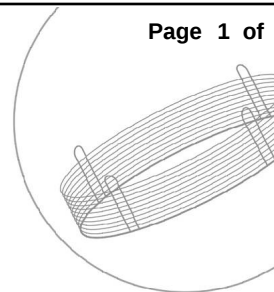

Products used in this Application:

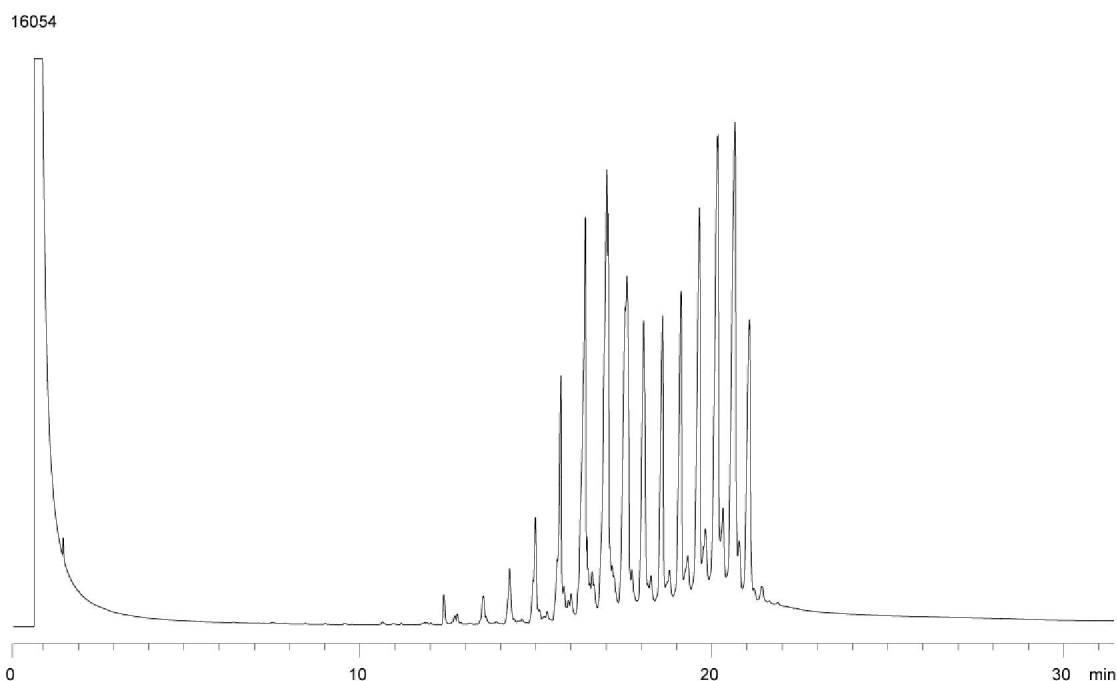Zebron™**ANALYTES:**

1 Butter

www.phenomenex.comPhenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.