

## Amino Acids

**Column:** Luna 3u C18(2) 100A

**Dimensions:** 150 x 2 mm ID

**Order No:** 00F-4251-B0

**Elution Type:** Isocratic


**Eluent A:** 0.1 M Sodium Acetate, pH 3.8 w/Acetic Acid

**Pct. A:** 100

**Flow Rate:** 0.12 mL/min

**Col. Temp.:** 40 °C


**Detection:** Fluorescence @ 385 nm / 515 nm (ambient)

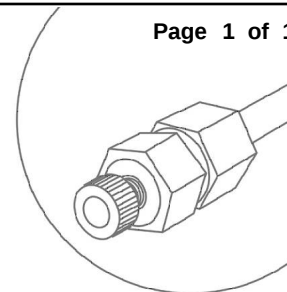

Page 1 of 1

Products used in this Application:

**LUNA**<sup>®</sup>

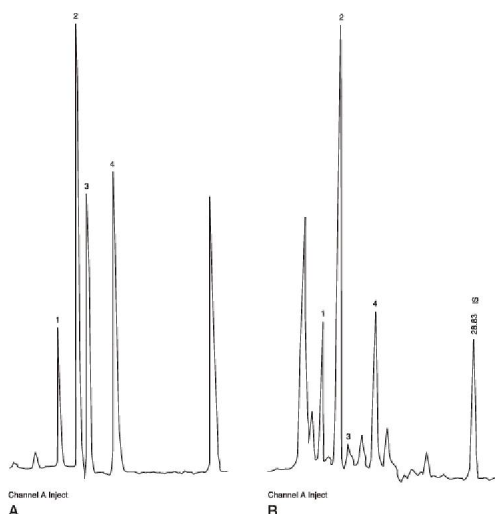

## ANALYTES:

- 1 Homocysteine
- 2 Cysteine
- 3 Cystinyl-glycine
- 4 Glutathione



[www.phenomenex.com](http://www.phenomenex.com)

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: [international@phenomenex.com](mailto:international@phenomenex.com).