

Glyphosate in Herbicide Roundup

Page 1 of 1

Column: Hamilton 10u PRP-X100 0A

Dimensions: 100 x 4.1 mm ID

Order No: CH0-0685

Elution Type: Isocratic



Eluent A: 40 mM Succinic acid pH 2.9

Pct. A: 100

Flow Rate: 1.7 mL/min

Col. Temp.: ambient

Detection: Conductivity (ambient)

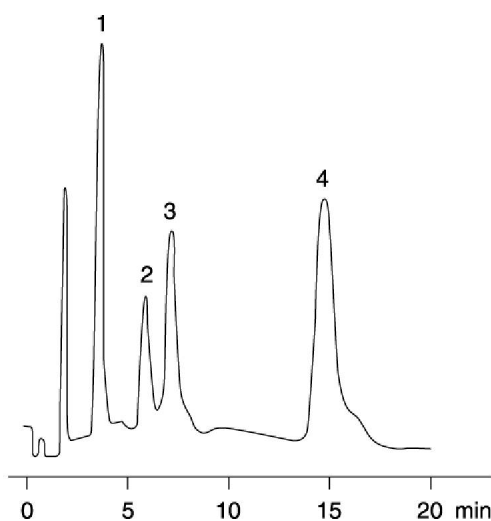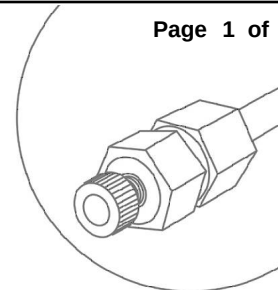

ANALYTES:

- 1 Fluoride
- 2 Glyphosate
- 3 Phosphate
- 4 Chloride



www.phenomenex.com

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.