

Antioxidants

Column: Ultracarb 5u ODS (30)

Dimensions: 150 x 4.6 mm ID

Order No: 00F-0351-E0

Elution Type: Gradient

Eluent A: Acetonitrile



Eluent B: Water

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	70	30
	2	10	100	0

Flow Rate: 1.7 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs. - Variable Wave. @ 280 nm (ambient)

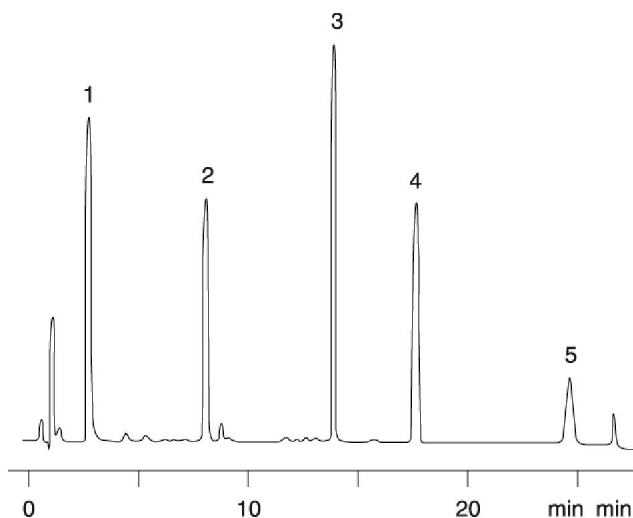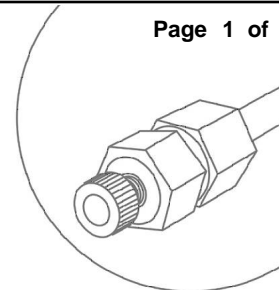

ANALYTES:

- 1 Ethyl 754
- 2 BHT
- 3 Goodrite 3125
- 4 Irganox 1010
- 5 Irganox 1076

