

Ginkgo Terpene Lactones and Flavone Glycosides Using Synergi Fusion-RP

Page 1 of 1

Column: Synergi 4u Fusion-RP 80A**Dimensions:** 150 x 4.6 mm ID**Order No:** 00F-4424-E0**Elution Type:** Gradient**Eluent A:** 0.1% Formic acid in Water**Eluent B:** 0.1% Formic acid in Acetonitrile


Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	100	0
	1	0	70	30
	2	15	30	70

Flow Rate: 1 mL/min**Col. Temp.:** ambient**Detection:** Mass Spectrometer (MS) (ambient)**Comment:** ESI positive

Products used in this Application:

synergi[™]

15737

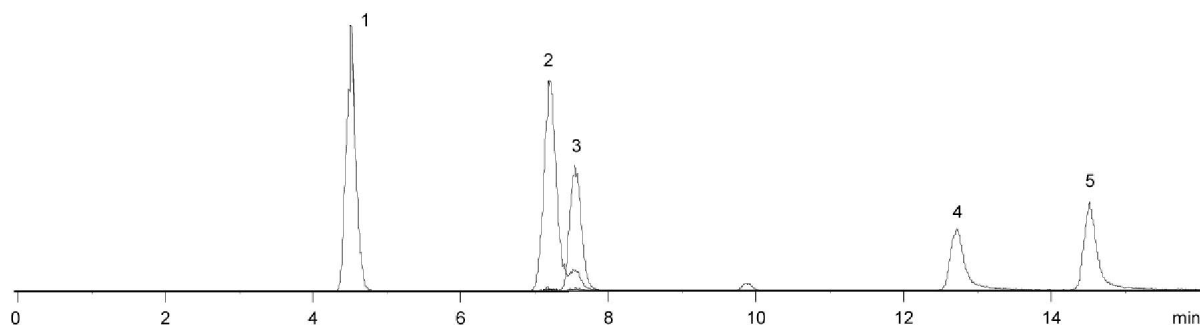

ANALYTES:

- 1 Bilobalide
- 2 Ginkgolide A
- 3 Ginkgolide B
- 4 Quercetin
- 5 Kaempferol



www.phenomenex.com

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.