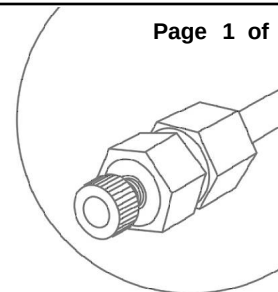
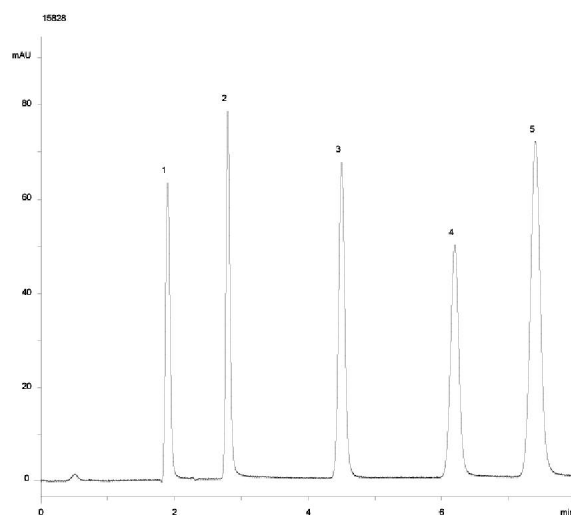


Nucleic Acid Bases on Gemini 3um C18 (W. Org Modifier)

Column: Gemini 3u C18 110A**Dimensions:** 100 x 4.6 mm ID**Order No:** 00D-4439-E0**Elution Type:** Isocratic**Eluent A:** 20mM Ammonium acetate, pH 4.5**Pct. A:** 99**Eluent B:** Acetonitrile**Pct. B:** 1**Flow Rate:** 1 mL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs. - Variable Wave. @ 260 nm (ambient)

Page 1 of 1

**Products used in this Application:****ANALYTES:**

- 1 Uracil
- 2 Cytosine
- 3 Guanine
- 4 Thymine
- 5 Adenine

www.phenomenex.comPhenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.