

Dexamethasone USP Method on Gemini 3um C18

Page 1 of 1

Column: Gemini 3u C18 110A

Dimensions: 100 x 4.6 mm ID

Order No: 00D-4439-E0

Elution Type: Isocratic

Eluent A: Water

Pct. A: 70


Eluent B: Acetonitrile

Pct. B: 30

Flow Rate: 1.5 mL/min

Col. Temp.: ambient


Detection: UV-Vis Abs. - Variable Wave. @ 254 nm (ambient)

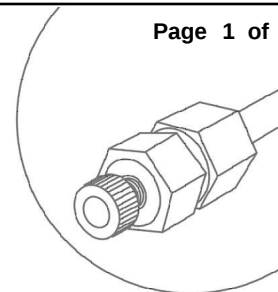

Products used in this Application:



Analyst Note: DIMENSIONS: 100 x 4.6mm, MOBILE PHASE: water:acetonitrile (70:30), TEMPERATURE: ambient, FLOW RATE: 1.5mL/min. DETECTION: 254 nm, SAMPLE: dexamethasone

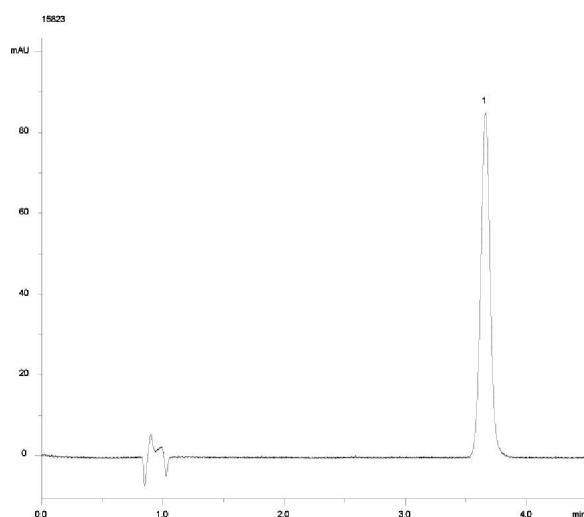

ANALYTES:

1 Dexamethasone

©Copyright 2000-2005 Phenomenex Inc.



www.phenomenex.com

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.