

Melamine and Impurities

Page 1 of 1

Column: Aqua 5u C18 125A

Dimensions: 250 x 4.6 mm ID

Order No: 00G-4299-E0

Elution Type: Isocratic


Eluent A: 20mM Potassium phosphate, pH 7

Pct. A: 100


Flow Rate: 1.5 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs. - Variable Wave. @ 254 nm (ambient)

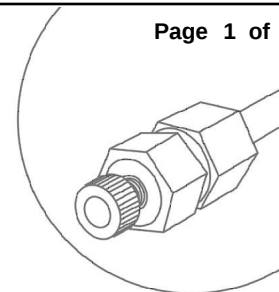

Products used in this Application:

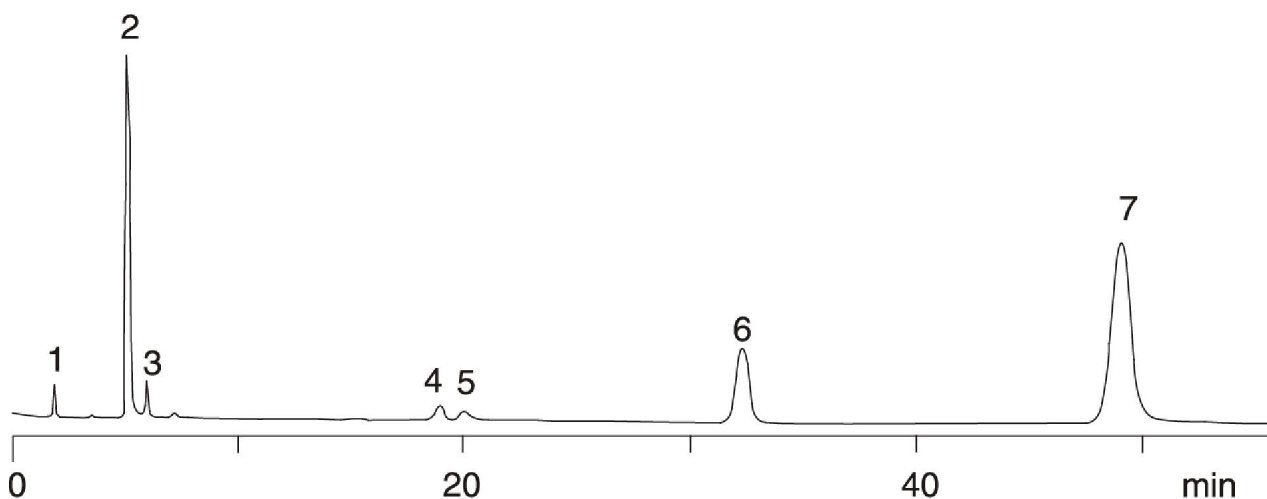AQUA[™]

ANALYTES:

- 1 Ammelide
- 2 Melamine
- 3 Cyanuric acid
- 4 Guanyl
- 5 Ureidomelamine
- 6 Melem
- 7 Melam