

Steroids

Column: Luna 5u C18(2) 100A

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4252-E0

Elution Type: Isocratic

Eluent A: 0.1% H₃PO₄ in Water

Eluent B: Acetonitrile

Eluent C: Methanol

Flow Rate: 1 mL/min

Col. Temp.: ambient


Detection: UV-Vis Abs. - Variable Wave. @ 254 nm (ambient)

Pct. A: 54

Pct. B: 35

Pct. C: 11


Page 1 of 1

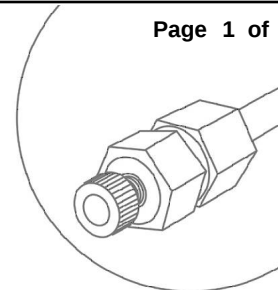

Products used in this Application:

LUNA®

14350

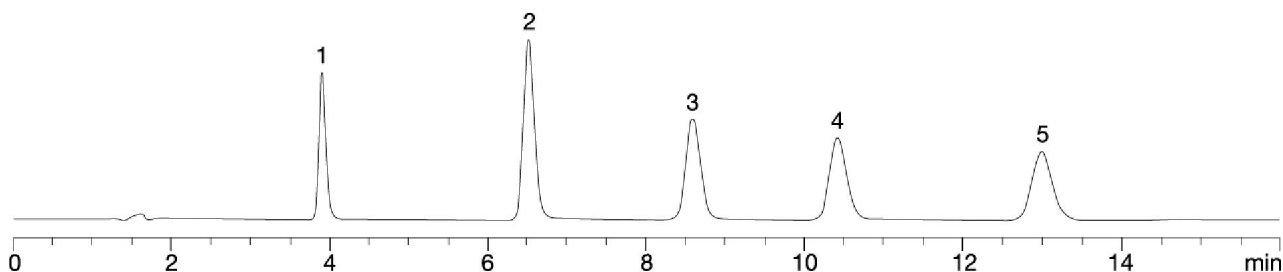

ANALYTES:

- 1 Hydrocortisone
- 2 Corticosterone
- 3 11- α -Hydroxy-Progesterone
- 4 Cortisone Acetate
- 5 11-Ketoprogesterone