


EPA Method 8081: Chlorinated Pesticides

Page 1 of 1


Column: Zebron ZB-50
Phase: 50%-diphenyl-50%-dimethyl polysiloxane copolymer
Dimensions: 30 meters x 0.25 mm x 0.25 µm
Order No: 7HG-G004-11
Oven Profile: 100°C for 0.2min to 200°C at 25°C/min to 240°C at 4°C/min to 320°C at 25°C/min for 2.5 min
Carrier Gas: Constant Flow Helium, 1.5 mL/min
Injection: Splitless 1 µL @ 250 °C
Detection: Electron Capture (ECD) (330 °C)

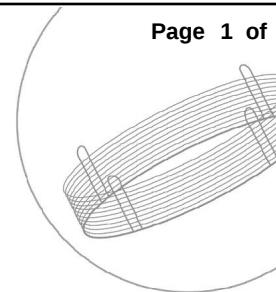

Products used in this Application:

Zebron[™]

15698

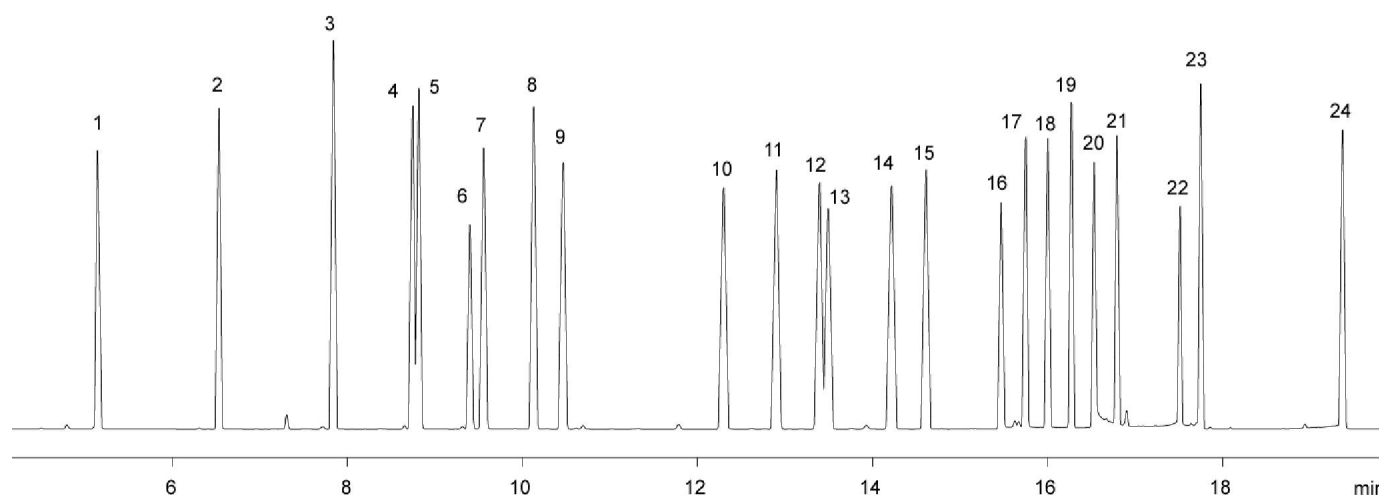

ANALYTES:

- | | | |
|----------------------------------|-----------------------|----------------------------|
| 1 1-Bromo-2-Nitrobenzene (IS) | 10 Heptachlor Epoxide | 19 Endrin Aldehyde |
| 2 Tetrachlorometaxylene (Surr) | 11 g-Chlordane | 20 Endosulfan Sulfate |
| 3 a-BHC | 12 a-Chlordane | 21 DDT |
| 4 b-BHC | 13 Endosulfan I | 22 Endrin Ketone |
| 5 g-BHC (Lindane) | 14 DDE | 23 Methoxychlor |
| 6 d-BHC | 15 Dieldrin | 24 Decachlorobiphenyl (IS) |
| 7 Pentachloronitrobenzene (Surr) | 16 Endrin | |
| 8 Heptachlor | 17 Endosulfan II | |
| 9 Aldrin | 18 DDD | |

©Copyright 2000-2005 Phenomenex Inc.

For more information call your Phenomenex Representative at (310) 212-0555



www.phenomenex.com

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.

phenomenex[®]
...breaking with tradition[™]

USA
tel.: (310) 212-0555
fax: (310) 328-7768
email: info@phenomenex.com

Puerto Rico
(800) 541-HPLC
(310) 328-7768
info@phenomenex.com

Canada
(800) 543-3681
(310) 328-7768
info@phenomenex.com

France
01 30 09 21 10
01 30 09 21 11
franceinfo@phenomenex.com

United Kingdom
01625-501367
01625-501796
ukinfo@phenomenex.com

Germany
06021-58830-0
06021-58830-11
anfrage@phenomenex.com

New Zealand
09-4780951
09-4780952
info@phenomenex.co.nz

Australia
02-9428-6444
02-9428-6445
info@phenomenex.com.au