

Capsaicin (Synergi 4µm Fusion-RP)

Column: Synergi 4µ Fusion-RP 80A

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4424-E0

Elution Type: Gradient

Eluent A: Water


Eluent B: Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	65	35
	2	29	40	60


Flow Rate: 1.5 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs. - Variable Wave. @ 227 nm (ambient)

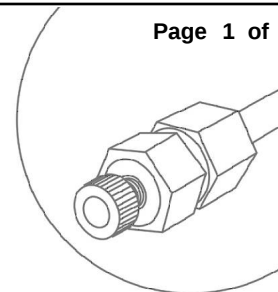

Products used in this Application:

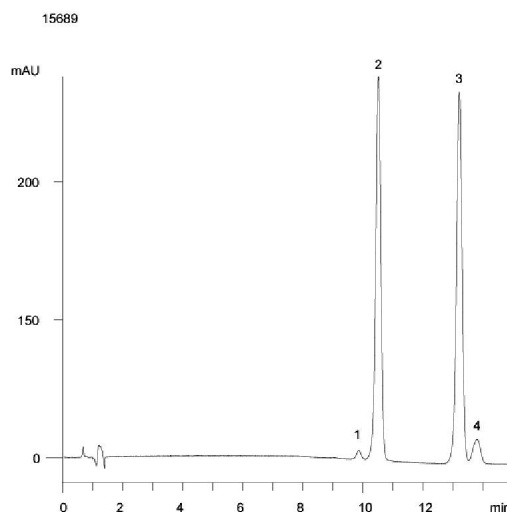synergi[™]

ANALYTES:

- 1 Nordihydrocapsaicin
- 2 Capsaicin
- 3 Dihydrocapsaicin
- 4 Homocapsaicin



www.phenomenex.com

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.