


Warfarin


Column: Chirex 3022 (S)-indoline-2-carboxylic acid and (R)-1-(alpha-naphthyl)ethylamine
Dimensions: 250 x 4.6 mm ID
Order No: 00G-3022-E0
Elution Type: Isocratic
Eluent A: Hexane/Dichloroethane/Ethanol-TFA (88:10:2)
Flow Rate: 1 mL/min
Col. Temp.: ambient
Detection: UV-Vis Abs. - Variable Wave. @ 282 nm (ambient)

Pct. A: 100

Page 1 of 1

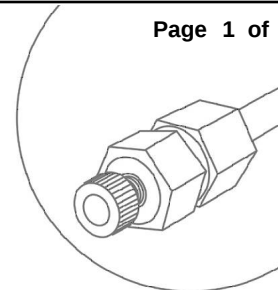

Products used in this Application:

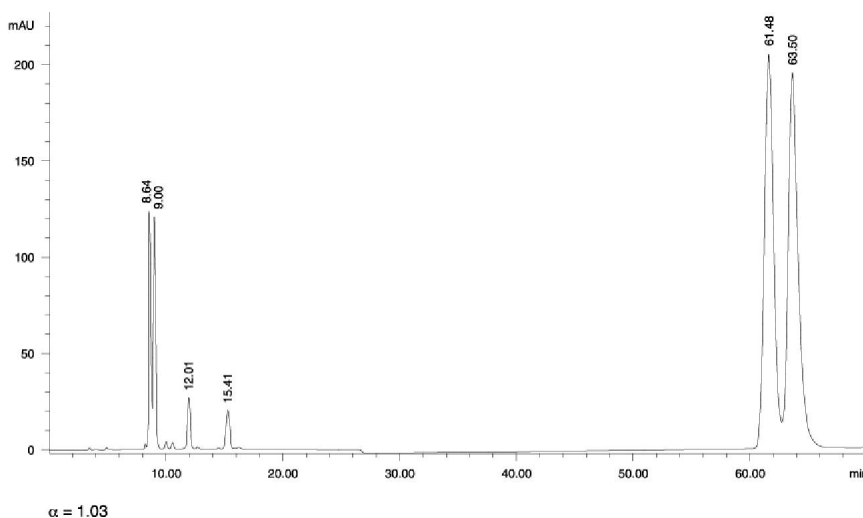CHIREX®

ANALYTES:

- 1 Warfarin
- 2 Warfarin