

Aromatic Polar Acids at pH 2.5 using Onyx C18 3.0mm ID

Page 1 of 1

Column: Onyx Monolithic C18

Dimensions: 100 x 3 mm ID

Order No: CH0-8158

Elution Type: Isocratic

Eluent A: 20mM Potassium phosphate pH 2.5

Pct. A: 75


Eluent B: Acetonitrile

Pct. B: 25

Flow Rate: 0.6 mL/min


Col. Temp.: ambient

Detection: UV-Vis Abs. - Variable Wave. @ 254 nm (ambient)

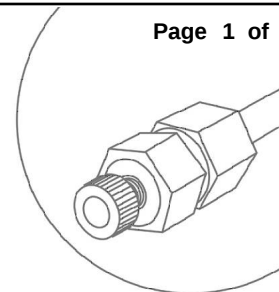

Products used in this Application:

ONYX

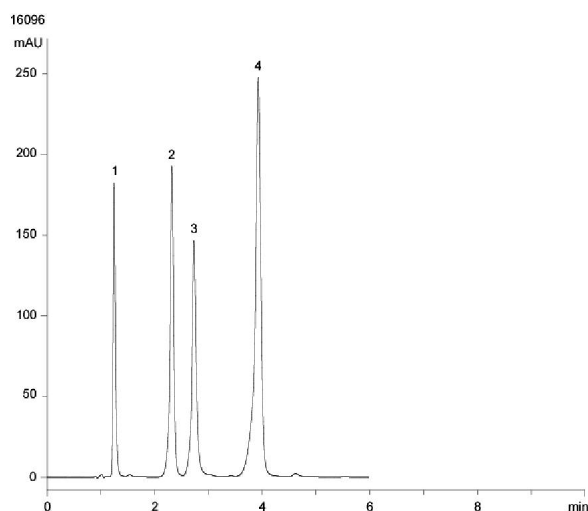

ANALYTES:

- 1 p-Hydroxybenzoic acid
- 2 Benzoic acid
- 3 Salicylic acid
- 4 p-Toluic acid