

Explosives

Column: Luna 5u Phenyl-Hexyl 100A

Dimensions: 250 x 4.6 mm ID

Order No: 00G-4257-E0

Elution Type: Isocratic

Eluent A: Methanol

Pct. A: 60

Eluent B: Water


Pct. B: 40

Flow Rate: 1.1 mL/min


Col. Temp.: ambient

Detection: UV-Vis Abs. - Variable Wave. @ 254 nm (ambient)

Page 1 of 1

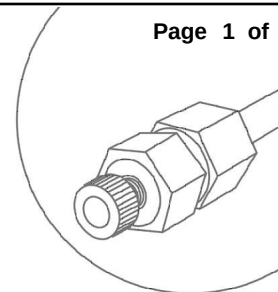

Products used in this Application:

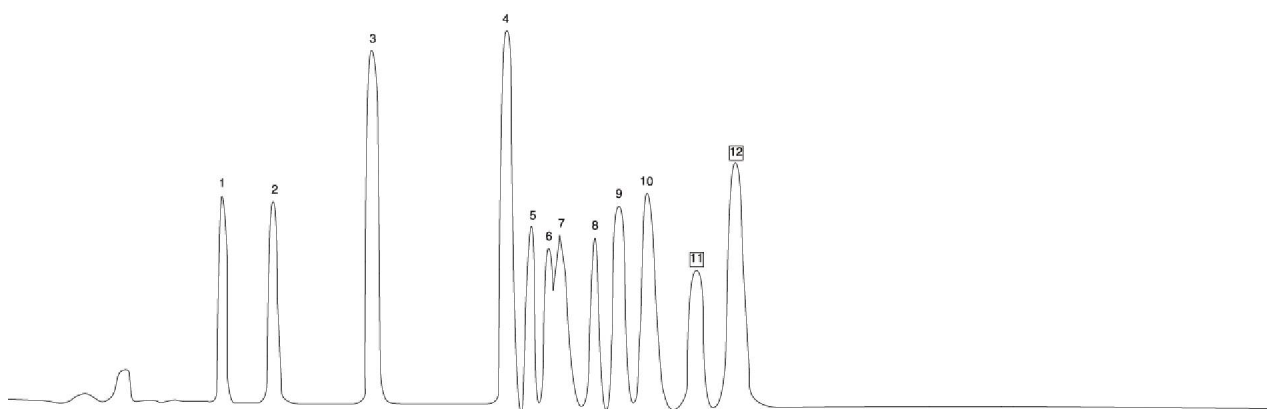LUNA®


ANALYTES:

- 1 HMX
- 2 RDX
- 3 NB
- 4 1,3-DNB
- 5 2-NT
- 6 4-NT
- 7 4-AM-DNT
- 8 3-NT
- 9 2-AM-DNT
- 10 1,3,5-TNB
- 11 2,6-DNT
- 12 2,4-DNT