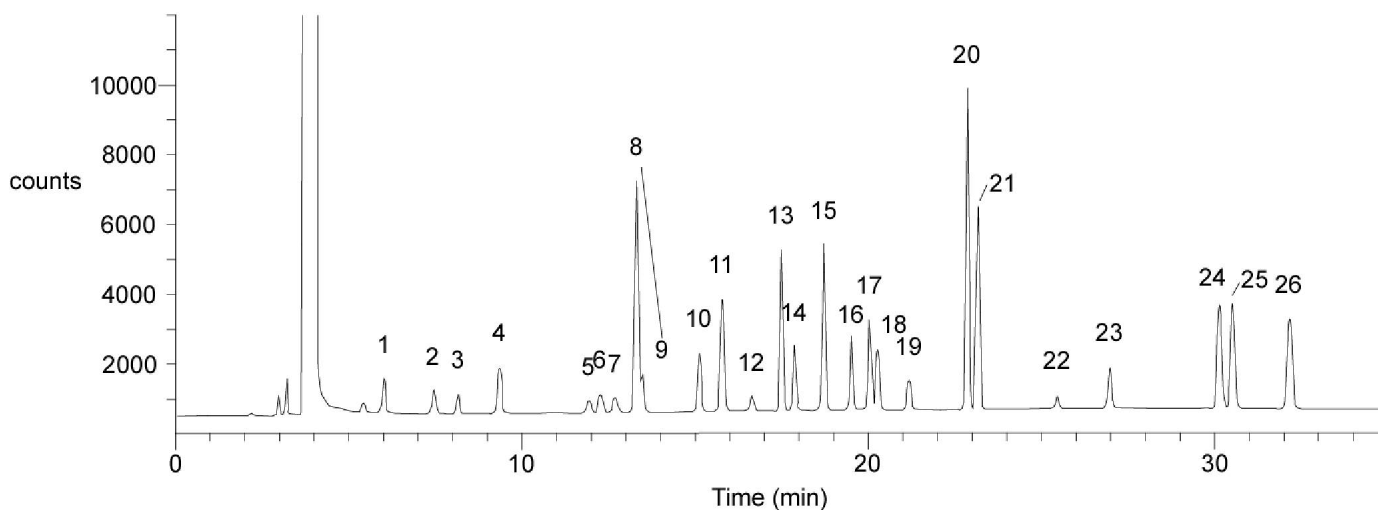


EPA Method 624 CLP

Page 1 of 1

Column: Zebron ZB-624
Phase: 6%-cyanopropylphenyl-94%-methyl polysiloxane
Dimensions: 30 meters x 0.53 mm x 3 µm
Order No: 7HK-G005-36
Oven Profile: 35°C (6 min), 5°C/min to 140°C (8 min)
Carrier Gas: Constant Flow Helium, 2.65 mL/min
Injection: Split 10:1 1 µL @ 250 °C
Detection: Flame Ionization (FID) (350 °C)

Products used in this Application:

ANALYTES:

- | | | |
|----------------------------|------------------------------|------------------------------|
| 1 1,1-Dichloroethene | 10 Trichloroethene | 19 Dibromochloromethane |
| 2 Methylene Chloride | 11 1,2-Dichloropropane | 20 Ethylbenzene |
| 3 trans-1,2-Dichloroethene | 12 Bromodichloromethane | 21 Chlorobenzene |
| 4 1,1-Dichloroethane | 13 2-Chloroethyl Vinyl Ether | 22 Bromoform |
| 5 Chloroform | 14 cis-1,3-Dichloropropene | 23 1,1,2,2-Tetrachloroethane |
| 6 1,1,1-Trichloroethane | 15 Toluene | 24 1,3-Dichlorobenzene |
| 7 Carbon Tetrachloride | 16 trans-1,3-Dichloropropene | 25 1,4-Dichlorobenzene |
| 8 Benzene | 17 1,1,2-Trichloroethane | 26 1,2-Dichlorobenzene |
| 9 1,2-Dichloroethane | 18 Tetrachloroethane | |

©Copyright 2000-2007 Phenomenex Inc.

For more information call your Phenomenex Representative at (310) 212-0555

www.phenomenex.com Phenomenex products are available worldwide. Email us at international@phenomenex.com					
Australia tel.: 02-9428-6444 email: info@phenomenex.com.au	Austria 01-319-1301 anfrage@phenomenex.com	Canada (800) 543-3681 info@phenomenex.com	Denmark 4824 8048 dkinfo@phenomenex.com	France 01 30 09 21 10 franceinfo@phenomenex.com	Germany 06021-58830-0 anfrage@phenomenex.com
Ireland tel.: 01 247 5405 email: eireinfo@phenomenex.com	Italy 051 736176 italiainfo@phenomenex.com	New Zealand 09-4780951 info@phenomenex.co.nz	Puerto Rico (800) 541-HPLC info@phenomenex.com	United Kingdom 01625-501367 ukinfo@phenomenex.com	USA (310) 212-0555 info@phenomenex.com