

Goldenseal Extract - Methanol Effects (cf. AppID 15471)

Page 1 of 1

Column: Luna 5u Phenyl-Hexyl 100A

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4257-E0

Elution Type: Isocratic

Eluent A: 20mM Potassium phosphate, pH 2.5

Pct. A: 50


Eluent B: Methanol

Pct. B: 50

Flow Rate: 1 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs. - Variable Wave. @ 254 nm (ambient)

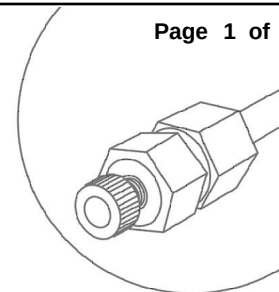

Products used in this Application:

LUNA[®]

15472



ANALYTES:

- 1 Hydrastatine
- 2 Berberine