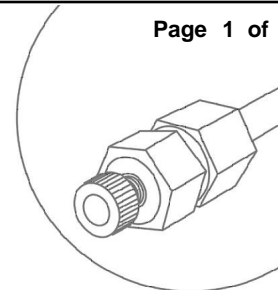
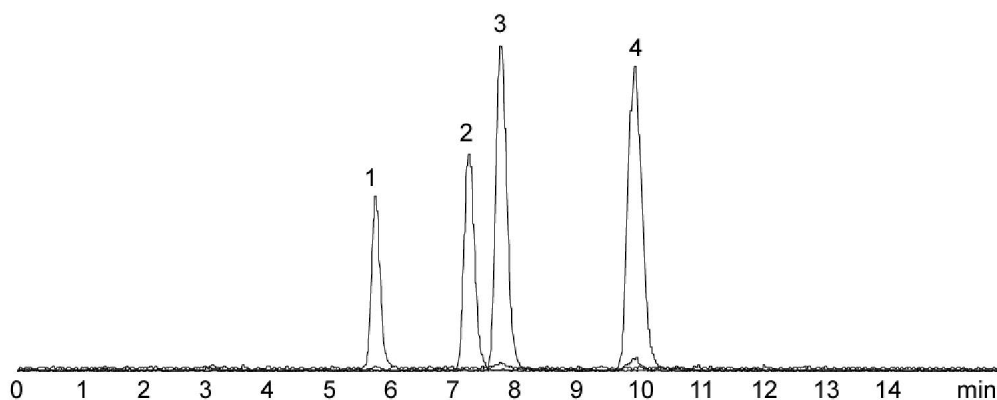


Aflatoxins using Synergi Max-RP

Column: Synergi 4 μ Max-RP 80Å
Dimensions: 150 x 3 mm ID
Order No: 00F-4337-Y0
Elution Type: Isocratic
Eluent A: 10 mM Ammonium Acetate, pH 4.5
Eluent B: Acetonitrile
Flow Rate: 0.5 mL/min
Col. Temp.: ambient
Detection: Mass Spectrometer (MS) (ambient)

Pct. A: 70**Pct. B:** 30**Products used in this Application:**synergi[™]

15523

**ANALYTES:**

- 1 Aflatoxin G2
- 2 Aflatoxin B2
- 3 Aflatoxin G1
- 4 Aflatoxin B1