

Citrinin analysis using Luna C18(2)

Page 1 of 1

Column: Luna 5u C18(2) 100A

Dimensions: 150 x 3 mm ID

Order No: 00F-4252-Y0

Elution Type: Isocratic

Eluent A: 0.1% Formic acid in Water

Pct. A: 50


Eluent B: Acetonitrile

Pct. B: 50

Flow Rate: 0.5 mL/min

Col. Temp.: ambient


Detection: Mass Spectrometer (MS) (ambient)

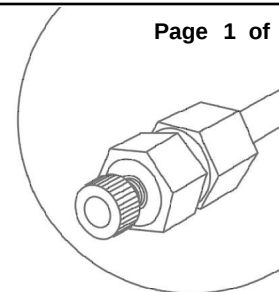

Products used in this Application:

LUNA®

15537

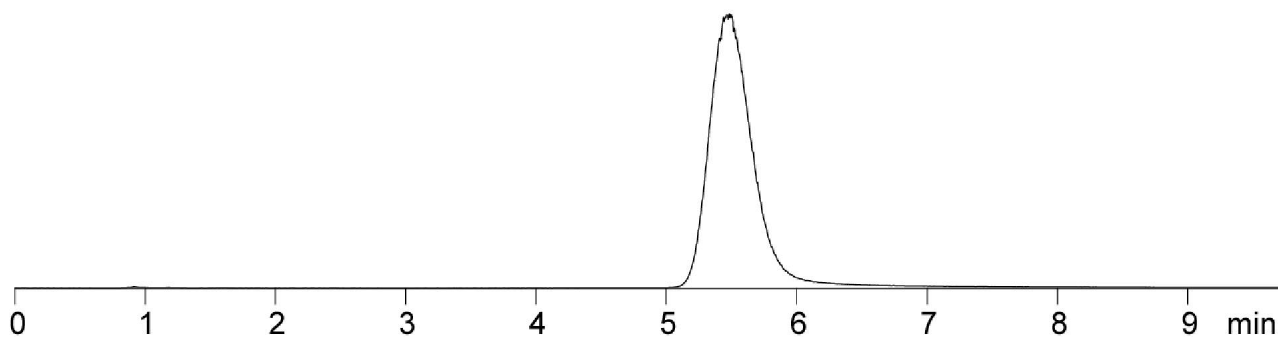

ANALYTES:

1 Citrinin