

Multivitamin tablet on Synergi Polar-RP

Page 1 of 1

Column: Synergi 4u Polar-RP 80A

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4336-E0

Elution Type: Gradient

Eluent A: 20mM KH₂PO₄ + 0.1% Hexane sulfonate, pH 3

Eluent B: Acetonitrile

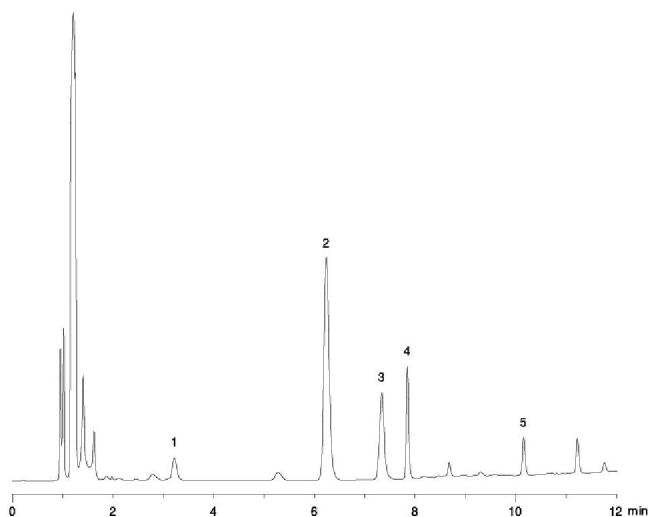
Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	97	3
	2	3	97	3
	3	18	50	50

Flow Rate: 1.5 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs. - Variable Wave. @ 210 nm (ambient)

Products used in this Application:

synergi[™]

ANALYTES:

- 1 Pantothenic acid
- 2 Pyridoxine
- 3 p-Aminobenzoic acid
- 4 Thiamine
- 5 Riboflavin



www.phenomenex.com

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.

USA
tel.: (310) 212-0555
fax: (310) 328-7768
email: info@phenomenex.com

Puerto Rico
(800) 541-HPLC
(310) 328-7768
info@phenomenex.com

Canada
(800) 543-3681
(310) 328-7768
info@phenomenex.com

France
01 30 09 21 10
01 30 09 21 11
franceinfo@phenomenex.com

United Kingdom
01625-501367
01625-501796
ukinfo@phenomenex.com

Germany
06021-58830-0
06021-58830-11
anfrage@phenomenex.com

New Zealand
09-4780951
09-4780952
info@phenomenex.co.nz

Australia
02-9428-6444
02-9428-6445
info@phenomenex.com.au