

Monosaccharides with organic acids using Rezex RHM

Page 1 of 1

Column: Rezex RHM-Monosaccharide H+(8%)

Dimensions: 300 x 7.8 mm ID

Order No: 00H-0132-K0

Elution Type: Isocratic


Eluent A: Water

Pct. A: 100

Flow Rate: 0.5 mL/min


Col. Temp.: 75 °C

Detection: Evaporative Light Scattering (ambient)

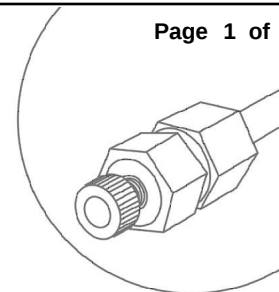

Products used in this Application:

REZEX™

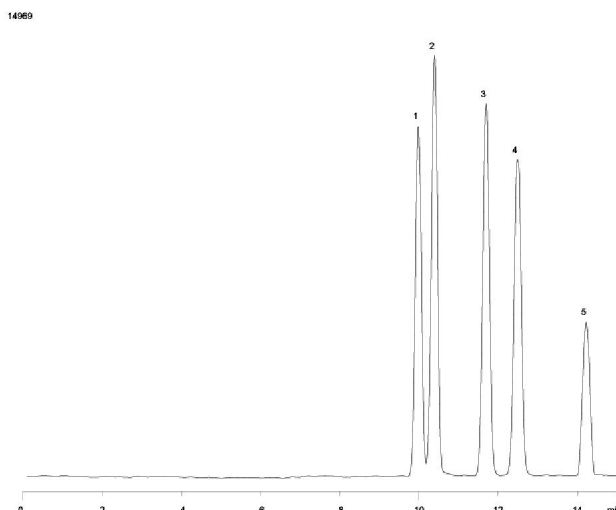

ANALYTES:

- 1 Citric acid
- 2 Tartaric acid
- 3 Glucose
- 4 Fructose
- 5 Succinic acid