



PONA: Retention Markers

Page 1 of 1

Column: Zebron ZB-1
Phase: 100% dimethyl polysiloxane
Dimensions: 100 meters x 0.25 mm x 0.5 µm
Order No: 7MG-G001-17
Oven Profile: 35°C (hold 14 min) to 60°C @ 1.1°C/min (hold 19 min) to 280°C @ 2°C/min (hold 5 min)
Carrier Gas: Constant Flow Helium, 1.55 mL/min
Injection: split 40:1 0.2 µL @ 220 °C
Detection: Flame Ionization (FID) (320 °C)

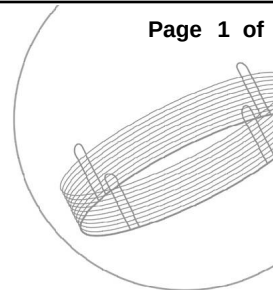

Products used in this Application:

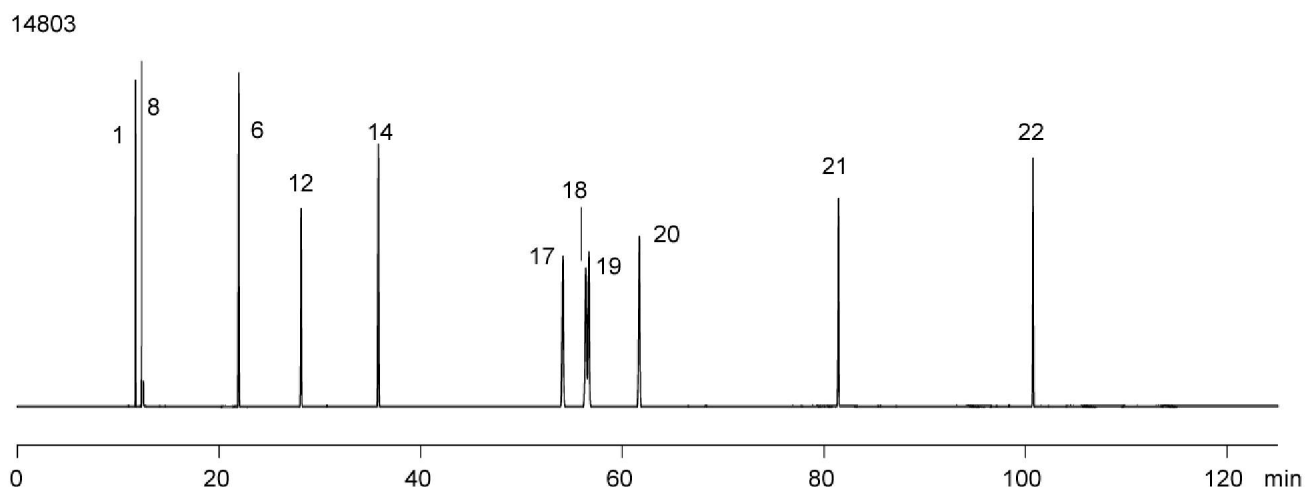Zebron™


ANALYTES:

- | | | |
|--------------------------|---------------------------|-------------|
| 1 Pentane | 10 3-Methylhexane | 19 p-Xylene |
| 2 2,2-Dimethylpentane | 11 2,2,4-Trimethylpentane | 20 o-Xylene |
| 3 Methylcyclopentane | 12 Heptane | 21 Decane |
| 4 Dimethylpentane | 13 2,3,4-Trimethylpentane | 22 Dodecane |
| 5 1-Methyl-1-Cyclohexane | 14 Toluene | |
| 6 Benzene | 15 2,2,3-Trimethylpentane | |
| 7 Cyclohexane | 16 2,3-Dimethylhexane | |
| 8 2-Methylhexane | 17 Ethylbenzene | |
| 9 2,3-Dimethylpentane | 18 m-Xylene | |

©Copyright 2000-2007 Phenomenex Inc.

For more information call your Phenomenex Representative at (310) 212-0555

www.phenomenex.com Phenomenex products are available worldwide. Email us at international@phenomenex.com					
Australia tel.: 02-9428-6444 email: info@phenomenex.com.au	Austria 01-319-1301 anfrage@phenomenex.com	Canada (800) 543-3681 info@phenomenex.com	Denmark 4824 8048 dkinfo@phenomenex.com	France 01 30 09 21 10 franceinfo@phenomenex.com	Germany 06021-58830-0 anfrage@phenomenex.com
Ireland tel.: 01 247 5405 email: eireinfo@phenomenex.com	Italy 051 736176 italiainfo@phenomenex.com	New Zealand 09-4780951 info@phenomenex.co.nz	Puerto Rico (800) 541-HPLC info@phenomenex.com	United Kingdom 01625-501367 ukinfo@phenomenex.com	USA (310) 212-0555 info@phenomenex.com