

Zearalone analysis using Gemini

Page 1 of 1

Column: Gemini 5u C18 110A
Dimensions: 150 x 3 mm ID
Order No: 00F-4435-Y0
Elution Type: Isocratic
Eluent A: 0.1% Formic acid in Water
Eluent B: Acetonitrile
Flow Rate: 0.5 mL/min
Col. Temp.: ambient
Detection: Mass Spectrometer (MS) (ambient)


Pct. A: 50

Pct. B: 50

Products used in this Application:

Gemini[™]

15533

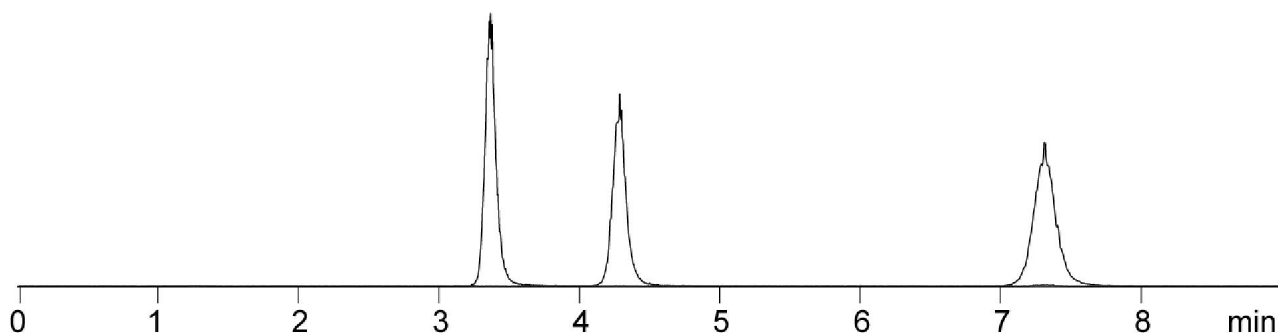

ANALYTES:

- 1 beta-Zearalenol
- 2 alpha-Zearalenol
- 3 Zearalenone



www.phenomenex.com

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.

USA
tel.: (310) 212-0555
fax: (310) 328-7768
email: info@phenomenex.com

Puerto Rico
(800) 541-HPLC
(310) 328-7768
info@phenomenex.com

Canada
(800) 543-3681
(310) 328-7768
info@phenomenex.com

France
01 30 09 21 10
01 30 09 21 11
franceinfo@phenomenex.com

United Kingdom
01625-501367
01625-501796
ukinfo@phenomenex.com

Germany
06021-58830-0
06021-58830-11
anfrage@phenomenex.com

New Zealand
09-4780951
09-4780952
info@phenomenex.co.nz

Australia
02-9428-6444
02-9428-6445
info@phenomenex.com.au