

Polar Acids

Column: Gemini 5u C18 110A

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4435-E0

Elution Type: Isocratic

Eluent A: 20mM Potassium Phosphat, pH 6.0

Pct. A: 70

Eluent B: Acetonitrile


Pct. B: 30

Flow Rate: 1 mL/min


Col. Temp.: 25 °C

Detection: UV-Vis Abs. - Variable Wave. @ 202 nm (25 °C)

Page 1 of 1

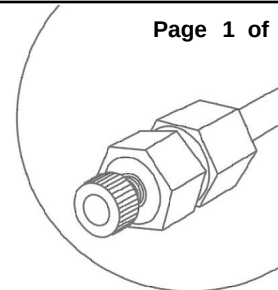

Products used in this Application:

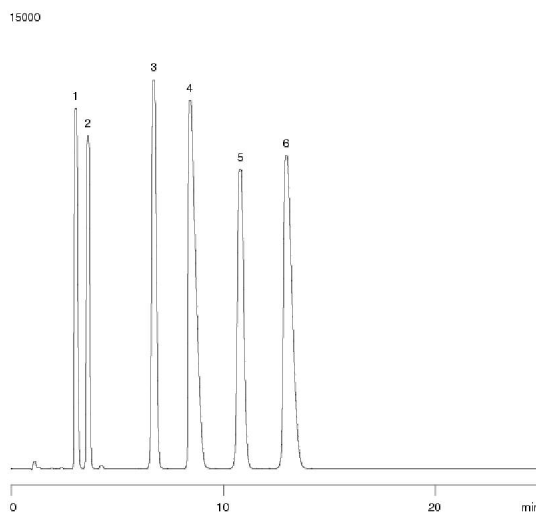

ANALYTES:

- 1 Tolmetin
- 2 Naproxen
- 3 Fenoprofen
- 4 Diflunisal
- 5 Indomethacin
- 6 Ibuprofen