

Cereal on Luna Amino

Column: Luna 5u NH2 100A

Dimensions: 250 x 4.6 mm ID

Order No: 00G-4378-E0

Elution Type: Isocratic

Eluent A: Acetonitrile

Eluent B: Water

Flow Rate: 40 mL/min


Col. Temp.: 40 °C

Detection: Refractive Index (RI) (ambient)


Pct. A: 80

Pct. B: 20

Page 1 of 1

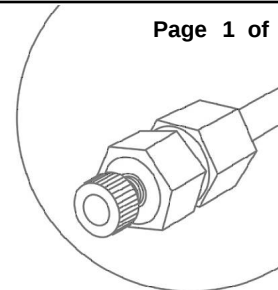

Products used in this Application:

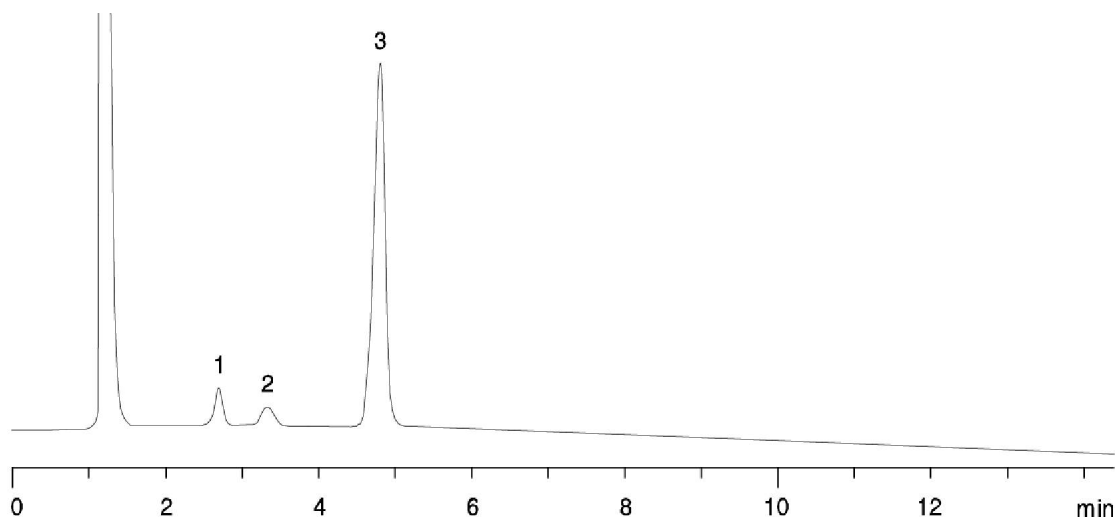LUNA®

ANALYTES:

- 1 Fructose
- 2 Glucose
- 3 Sucrose