

Apple Juice analysis using Rezex RCM

Column: Rezex RCM-Monosaccharide Ca+ (8%)

Dimensions: 300 x 7.8 mm ID

Order No: 00H-0130-K0

Elution Type: Isocratic

Eluent A: Water


Flow Rate: 0.6 mL/min

Col. Temp.: 80 °C

Detection: Evaporative Light Scattering (ambient)


Pct. A: 100

Page 1 of 1

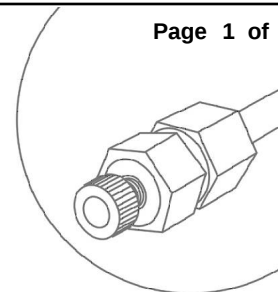

Products used in this Application:

REZEX™

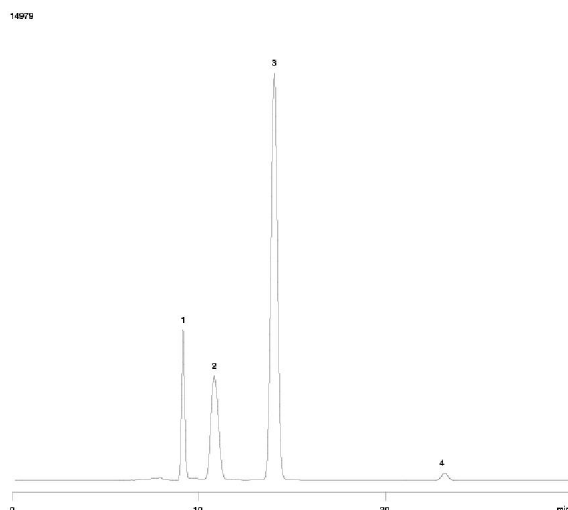

ANALYTES:

- 1 Sucrose
- 2 Glucose
- 3 Fructose
- 4 Sorbitol