

Polyethylene Glycol (PEG)

Column: Unknown or Proprietary Material

Dimensions: 300 x 8 mm ID

Order No:

Elution Type: Isocratic

Eluent A: Water



Pct. A: 100

Flow Rate: 1 mL/min

Col. Temp.: ambient

Detection: Refractive Index (RI) (ambient)

Page 1 of 1

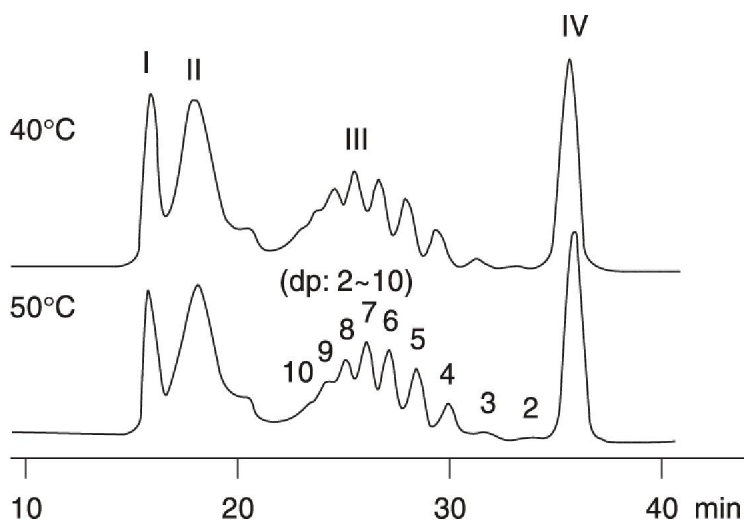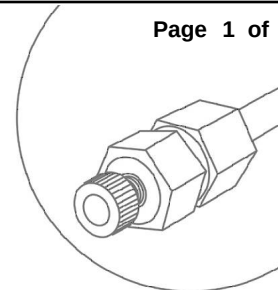

ANALYTES:

- 1 Polyethylene glycol (PEG) MW 4000
- 2 Polyethylene glycol (PEG) MW 1000
- 3 Polyethylene glycol (PEG) MW 300
- 4 Ethylene glycol (EG)