


Organic acids on Synergi Hydro-RP

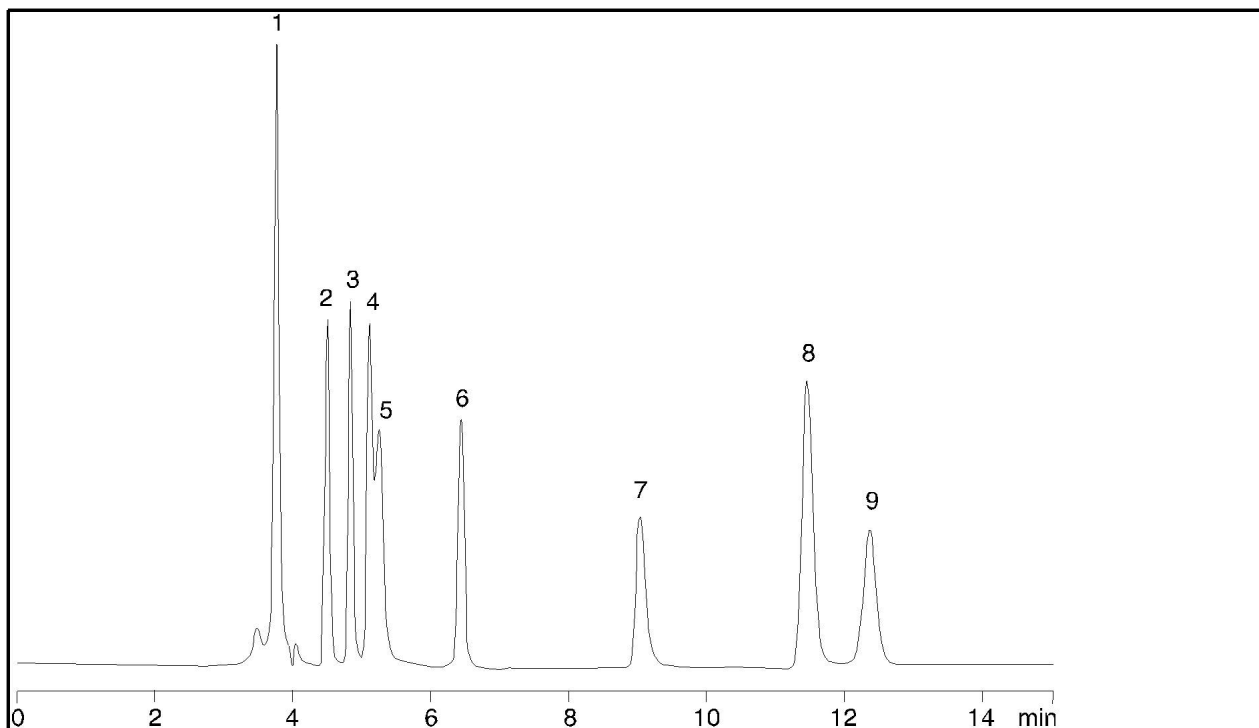

Column: Synergi 4 μ m C18 / ODS Hydro-RP

Dimensions: 250 x 4.6 mm ID

Order No:

Elution Type: Isocratic

Eluent A: 20mM Potassium phosphate pH 2.9

Pct. A: 100

Flow Rate: 0.7 mL/min

Temperature: 22 °C

Detection: UV-Vis Abs. - Variable Wave. @ 220 nm

SAMPLE ANALYTE

- 1 Oxalic acid
- 2 Tartaric acid
- 3 Glycolic acid
- 4 Formic acid
- 5 Pyruvic acid
- 6 Malonic acid
- 7 Acetic acid
- 8 Maleic acid
- 9 Citric acid

SPE:

For more information call your Phenomenex Representative at (310) 212-0555

Copyright 2000 Phenomenex Inc.

Ref. No.: 14270