


Tocopherols (Vitamin E) on Luna Amino

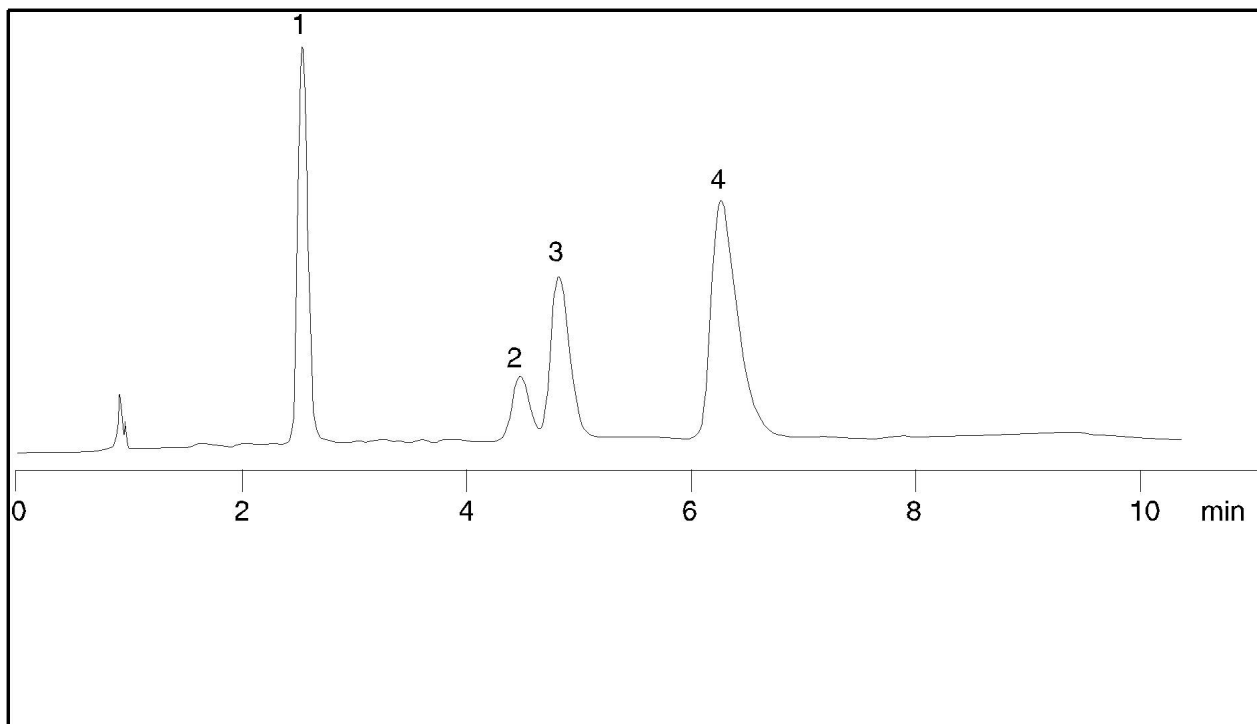

Column: Luna 5 μ m NH2

Dimensions: 250 x 4.6 mm ID

Order No:

Elution Type: Isocratic

Eluent A: Hexane

Pct. A: 97

Eluent B: Isopropyl Alcohol

Pct. B: 3

Flow Rate: 4 mL/min

Temperature: 22 °C

Detection: UV-Vis Abs. - Variable Wave. @ 298 nm

SAMPLE ANALYTE

- 1 *alpha*-Tocopherol
- 2 *beta*-Tocopherol
- 3 *gamma*-Tocopherol
- 4 *delta*-Tocopherol

SPE:

For more information call your Phenomenex Representative at (310) 212-0555

Copyright 2000 Phenomenex Inc.

Ref. No.: 14296