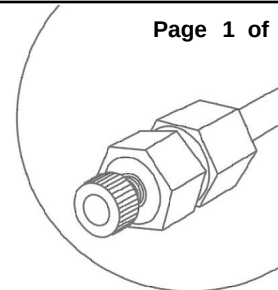
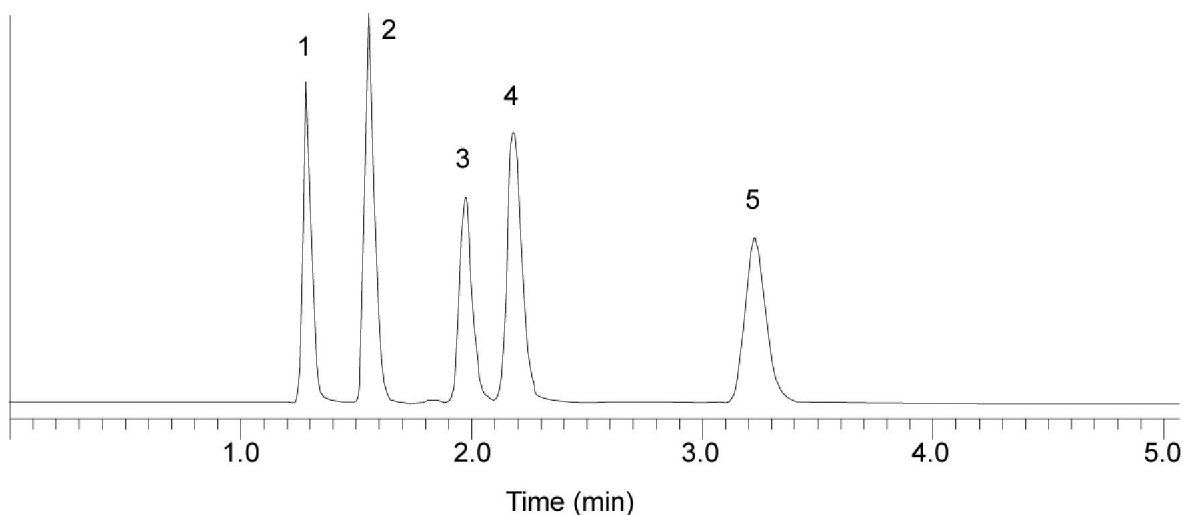


Nucleic Acid Bases

Column: Luna 5u NH2 100A**Dimensions:** 150 x 4.6 mm ID**Order No:** 00F-4378-E0**Elution Type:** Isocratic**Eluent A:** Acetonitrile**Pct. A:** 80**Eluent B:** Water**Pct. B:** 20**Flow Rate:** 1 mL/min**Col. Temp.:** 40 °C**Detection:** UV-Vis Abs. - Variable Wave. @ 254 nm (ambient)

Page 1 of 1

**Products used in this Application:****LUNA[™]****ANALYTES:**

- 1 Thymine
- 2 Uracil
- 3 Cytosine
- 4 Adenine
- 5 Guanosine

**www.phenomenex.com**

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.