

## Cytochrome C Bacteria

Column: Luna 3u C18(2) 100A

Dimensions: 150 x 0.5 mm ID

Order No: 00F-4251-AF

Elution Type: Gradient

Eluent A: 0.01% TFA in water


Eluent B: 0.008% TFA in CH<sub>3</sub>CN

| Gradient Profile: | Step No. | Time (min) | Pct A | Pct B |
|-------------------|----------|------------|-------|-------|
|                   | 1        | 0          | 95    | 5     |
|                   | 2        | 60         | 40    | 60    |
|                   | 3        | 60.01      | 95    | 5     |
|                   | 4        | 70         | 95    | 5     |


Flow Rate: 15 µL/min

Col. Temp.: 40 °C

Detection: UV-Vis Abs. - Variable Wave. @ 210 nm (ambient)

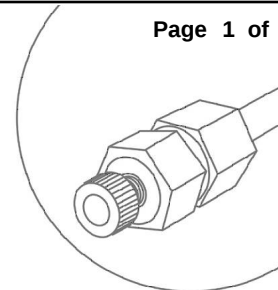

Products used in this Application:

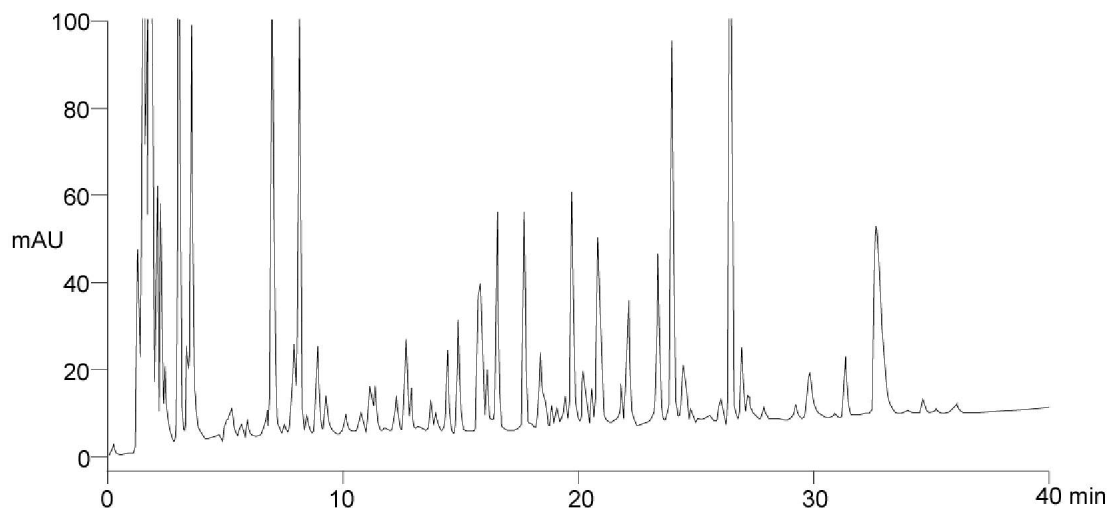**LUNA**<sup>™</sup>

## ANALYTES:

- 1 Cythochrom C bacteria tryptic digest



www.phenomenex.com

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.