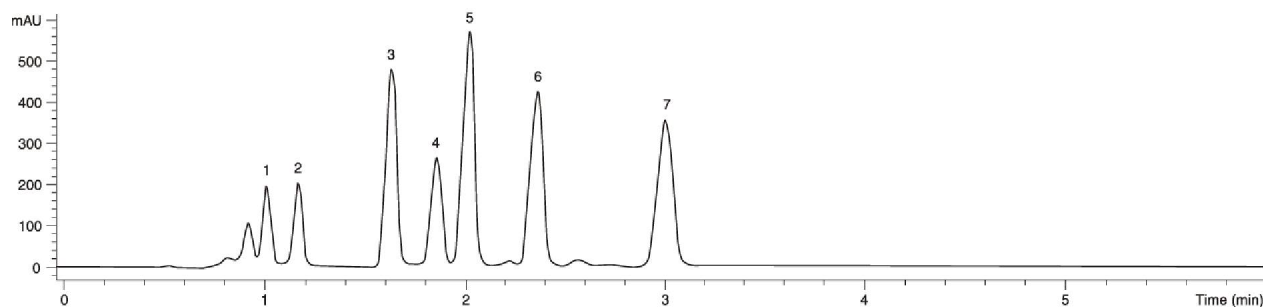


Quinones

LUNA™

G A L A X Y™
S E R I E S**Column:** Luna 3 µm C18 / ODS (2)**Dimensions:** 75 x 4.6 mm ID**Order No:** 00C-4251-E0**Elution Type:** Isocratic**Eluent A:** Methanol**Pct. A:** 89**Eluent B:** Water**Pct. B:** 11**Flow Rate:** 1 mL/min**Temperature:** 22 °C**Detection:** UV-Vis Abs. - Variable Wave. @ 251 nm**SAMPLE ANALYTES:**

- 1 1, 4-Naphthoquinone
- 2 2-Methyl-1, 4-naphthoquinone
- 3 9,10-Anthroquinone
- 4 1,5-Dichloro-9,10-anthaquinone
- 5 2-Methyl-9,10-anthroquinone
- 6 2-Ethyl-9,10-anthroquinone
- 7 2-tert-Butyl-9,10-anthroquinone

SPE:**For more information call your Phenomenex Representative at (310) 212-0555**

©Copyright 2000-2005 Phenomenex Inc.

Ref. No.: 1188/ALG/39

Page 1 of 1