

Tryptic Digest of RNase (deamidation) on Jupiter Proteo

Page 1 of 1

Column: Jupiter 4u Proteo 90A

Dimensions: 250 x 4.6 mm ID

Order No: 00G-4396-E0

Elution Type: Gradient

Eluent A: 1.25mM TFA in Water

Eluent B: 1.0mM TFA in Acetonitrile

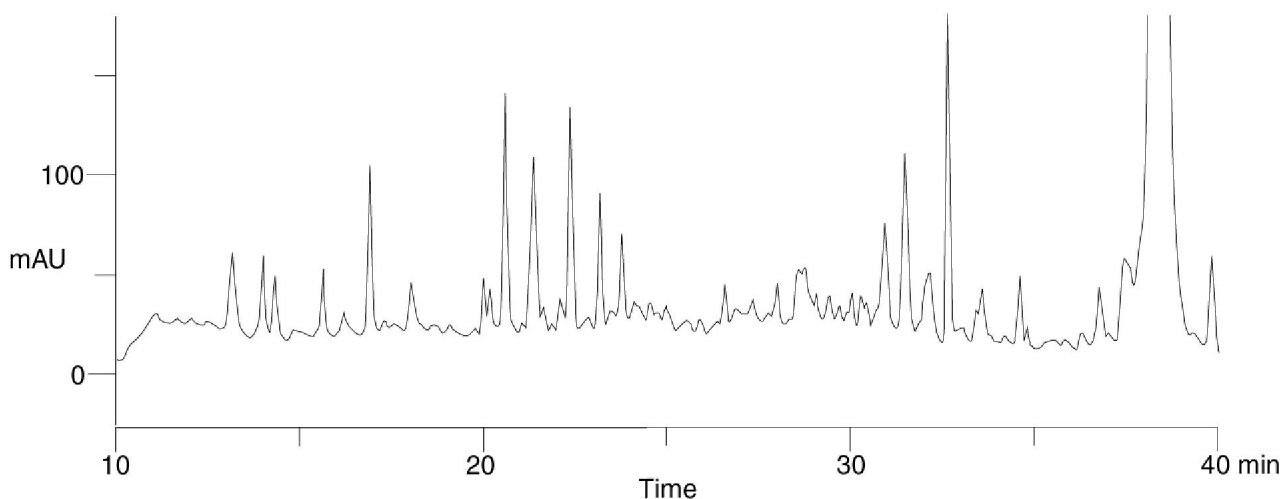
Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	95	5
	2	5	95	5
	3	60	60	40

Flow Rate: 1 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs. - Variable Wave. @ 210 nm (ambient)

Products used in this Application:

Jupiter

ANALYTES:

- 1 Tryptic Digest of RNase



www.phenomenex.com

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.