


Tricyclic Antidepressants at pH 11 in Isocratic RP HPLC



15016

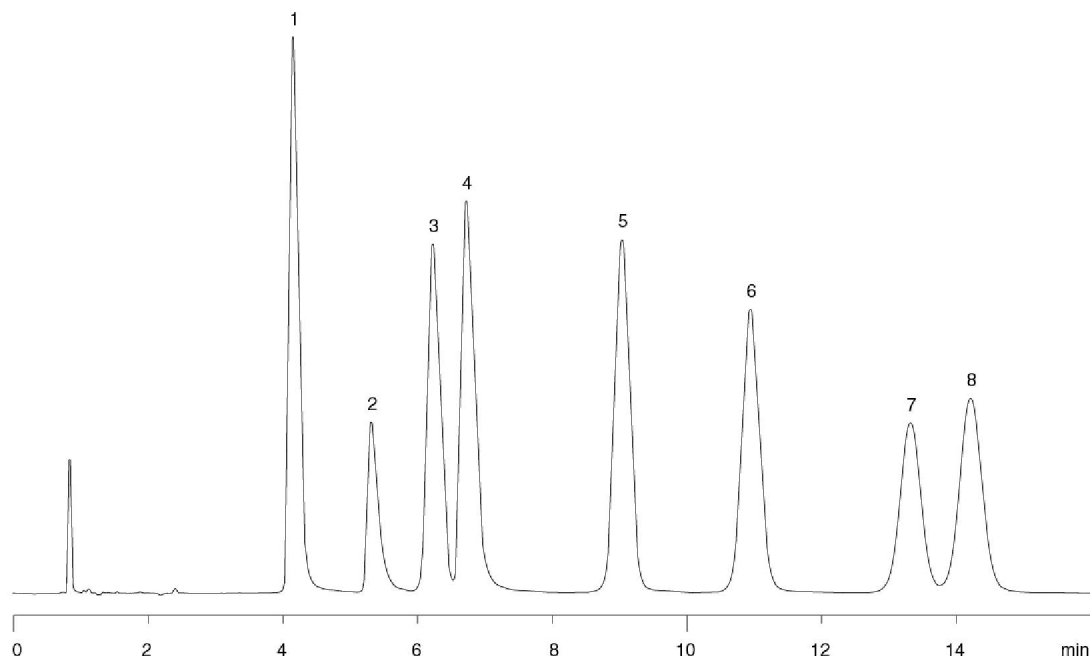

Column: Gemini 5 μ m C18 / ODS C18

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4435-E0

Elution Type: Isocratic

Eluent A: 20mM TEA, pH11.0

Pct. A: 30

Eluent B: ACN

Pct. B: 25

Eluent C: MeOH

Pct. C: 45

Flow Rate: 1.5 mL/min

Temperature: 22 °C

Detection: UV-Vis Abs. - Variable Wave. @ 254nm nm

SAMPLE ANALYTES:

- 1 Nordoxepine
- 2 Protriptyline
- 3 Doxepine
- 4 Nortriptyline
- 5 Imipramine
- 6 Amitriptyline
- 7 Clomipramine
- 8 Trimipramine

SPE:

For more information call your Phenomenex Representative at

©Copyright 2000-2005 Phenomenex Inc.

Ref. No.: 15016

Page 1 of 1