


APPLICATION NOTE - HPLC
Peptides on Jupiter Proteo

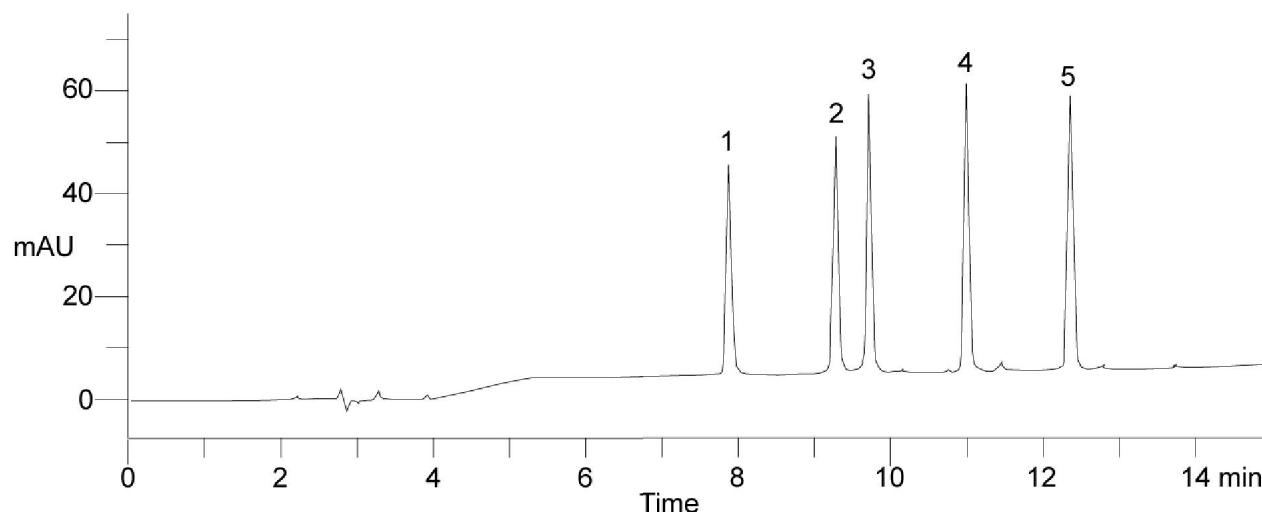

Column: Jupiter 4 μ m C12 Proteo
Dimensions: 250 x 4.6 mm ID
Order No: 00G-4396-E0
Elution Type: Gradient
Eluent A: 0.1% TFA in Water
Eluent B: 0.08% TFA in Acetonitrile

Gradient Profile:

Step No.	Step Time	Pct A	Pct B
1	0	95	5
2	20	60	40

Flow Rate: 1 mL/min
Temperature: 22 °C
Detection: UV-Vis Abs. - Variable Wave. @ 210 nm
SPE:

SAMPLE ANALYTES:

- 1 Arg-Gly-Ala-Gly-Gly-Leu-Gly-Leu-Gly-Lys-amide
- 2 Ac-Arg-Gly-Gly-Gly-Gly-Leu-Gly-Leu-Gly-Lys-amide
- 3 Ac-Arg-Gly-Ala-Gly-Gly-Leu-Gly-Leu-Gly-Lys-amide
- 4 Ac-Arg-Gly-Val-Gly-Gly-Leu-Gly-Leu-Gly-Lys-amide
- 5 Ac-Arg-Gly-Val-Val-Gly-Leu-Gly-Leu-Gly-Lys-amide

For more information call your Phenomenex Representative at (310) 212-0555

©Copyright 2000-2005 Phenomenex Inc.

Ref. No.: 14396

Page 1 of 1