


APPLICATION NOTE - HPLC
Polar Bases -- Beta-Blockers



15023

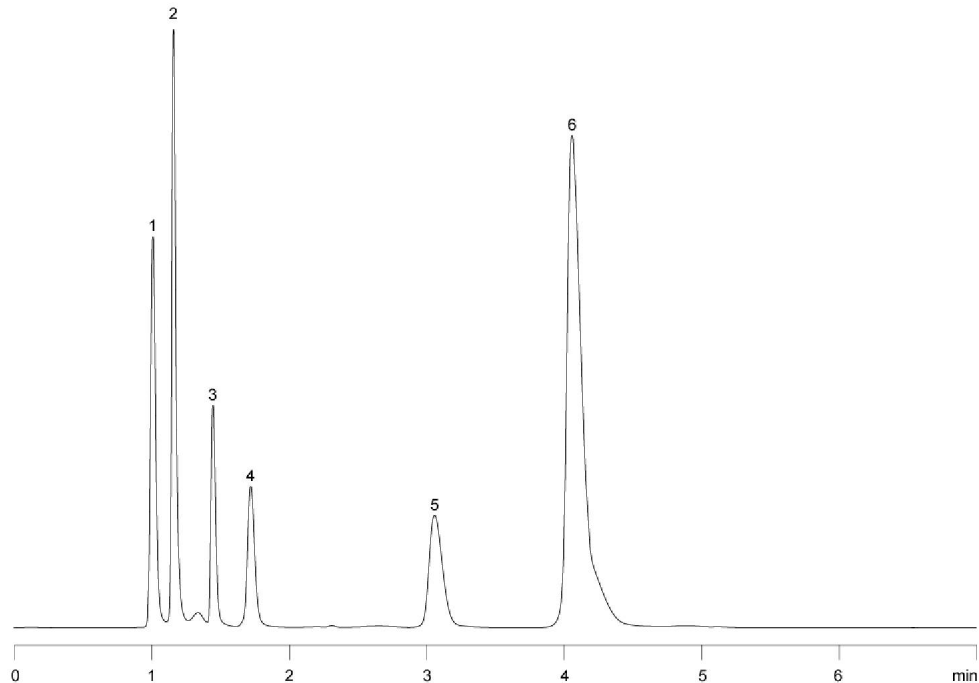

Column: Gemini 5 μ m C18 / ODS C18

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4435-E0

Elution Type: Isocratic

Eluent A: 20mM Phosphate Buffer, pH2.5

Pct. A: 75

Eluent B: ACN

Pct. B: 25

Flow Rate: 1.5 mL/min

Temperature: 22 $^{\circ}$ C

Detection: UV-Vis Abs. - Variable Wave. @ 230 nm

SAMPLE ANALYTES:

- 1 *Atenolol*
- 2 *Pindolol*
- 3 *Metoprolol*
- 4 *Labetalol*
- 5 *Propranolol*
- 6 *Alprenolol*

SPE:

For more information call your Phenomenex Representative at

©Copyright 2000-2005 Phenomenex Inc.

Ref. No.: 15023

Page 1 of 1

